

# Specialty Cabinets

Safe and Secure







**Fast & Professional**

## About our Specialty Cabinets

Our cabinets are designed for specialized storage needs with durability and security in mind. The Acid/Corrosive Resistant Vented Safety Cabinets provide robust storage for acids and solvents, featuring adjustable shelves and optional ADA-compliant height. Flammable cabinets come UL rated with manual or self-closing options, adjustable shelf and lock. Both cabinets offer customizable features like door swings and finished backs to meet specific facility requirements.

## Options

-  **Acid/Corrosive Storage Cabinet**
-  **Flammable Storage Cabinet**
-  **Narcotic Locker Cabinet**
-  **Vacuum Pump Storage**



[www.FENCOLABS.com](http://www.FENCOLABS.com)

Find your solution today!



Toll-Free: 800.486.8484

# SPECIALTY CABINETS


## ACID & FLAMMABLE STORAGE

### Product Features:

- Cabinets feature the same construction as standard cabinets with additional features as specified.
- Shelves are adjustable in 1/2" increments

### Ordering Information

- For finished backs, add suffix 'FB' to model number.
- ADA Height unit also available, 32-1/8"H

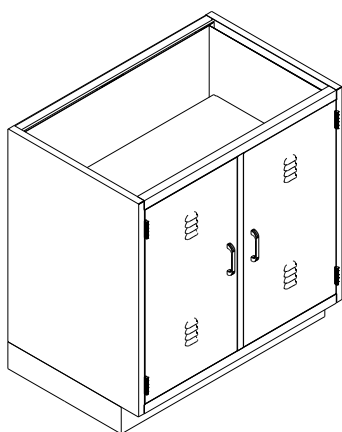

### Acid/Corrosive Resistant Vented Safety Cabinet


The 7000 cabinet safely stores acids within a welded steel unit lined with Wesliner material for maximum corrosion resistance in the event of spills or accidents. Each of the two hinged doors is vented top and bottom to prevent any build up of noxious fumes

The **Phenolic Resin Liner** used is 1/8" thick with a natural white smooth surface on both sides. It is a strong and rugged material designed to offer maximum resistance to heat, chemicals, impact, and moisture.

Unit is standard 35-3/4"H x 22-1/2"D and includes one full width white Phenolic resin adjustable shelf. Standard "ACID" label is also included (red letters on silver).

### Modifications for 7000:

- Specify swing for cupboard doors on **single** units by substituting '\*' with 'R' or 'L' after series number (Units are shown with left swing).  
Ex. 7009R24 for right swing cupboard unit, 24" Wide)

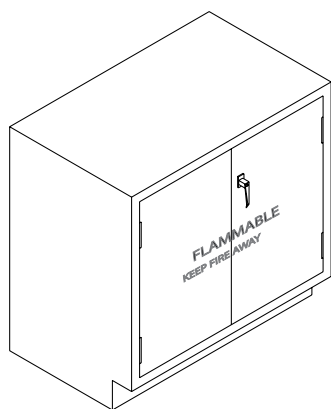

### Flammable Storage Cabinet

Cabinets for storing flammables in safety cans/containers.

#### Manual Close Doors

<b>CF3848</b> .....	24"W 35-3/4"H 22"D
<b>CF3847</b> .....	30"W 35-3/4"H 22"D
<b>CF3845</b> .....	35"W 35-3/4"H 22"D
<b>CF4000</b> .....	35"W 35-9/16"H 18"D
<b>L124</b> .....	35"W 35-9/16"H 22"D
<b>A145</b> .....	43"W 65"H 18"D
<b>CF3836</b> .....	47"W 35-3/4"H 22"D
<b>L136</b> .....	47"W 35-9/16"H 22"D
<b>CF3884</b> .....	48"W 35-3/4"H 18"D
<b>L144</b> .....	59"W 35-9/16"H 22"D
<b>CF3885</b> .....	35"W 35.75"H 18"D

#### Self Close Doors

<b>CF3882</b> .....	24"W 35-3/4"H 22"D
<b>CF3880</b> .....	30"W 35-3/4"H 22"D
<b>CF3881</b> .....	35"W 35-3/4"H 22"D
<b>L224</b> .....	35"W 35-5/8"H 22"D
<b>A245</b> .....	43"W X 65"H X 18"D
<b>CF3878</b> .....	47"W 35-3/4"H 22"D
<b>L236</b> .....	47"W X 35-5/8"H 22"D
<b>CF3883</b> .....	48"W 35-3/4"H 22"D
<b>L244</b> .....	59"W 35-5/8"H 22"D

### Acid/Corrosive Cabinet

Cabinets for storing of corrosives/flammables/acids.

#### Manual Close Doors

<b>C124</b> .....	35"W 36"H 22"D
<b>CC3844</b> .....	35"W 35-3/4"H 22"D
<b>C145</b> .....	43"W 65"H 18"D
<b>CF3886</b> .....	30"W 35.75"H 22"D

#### Self Close Doors

<b>CC3879</b> .....	35"W 35-3/4"H 22"D
<b>C245</b> .....	43"W 65"H 18"D
<b>CC3885</b> .....	47"W 35-3/4"H 22"D
<b>C224</b> .....	35"W 36"H 22"D

# SPECIALTY CABINETS


## NARCOTICS STORAGE CABINETS

**Product Features:**

- Sound deadening material
- Integral door and toe kick; no bottom of cabinet
- Cutout for electrical switch on front and rear
- Pullout lower shelf options available

**Ordering Information**

- Specify swing for cupboard doors on single units by substituting '\*' with 'R' or 'L' after series number (Units are shown with left swing).
- For pullout shelf, add suffix "-PO" to end of model number (7019L24VAC-PO)

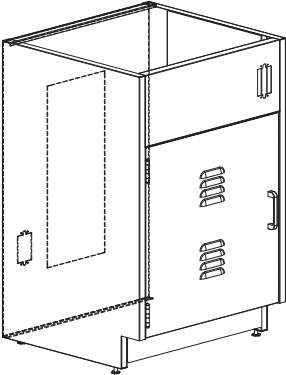


**Vacuum Pump Cabinets**

- Sound deadening material
- Integral door and toe kick; no bottom of cabinet
- Cutout for electrical switch on front and rear
- Pullout lower shelf options available

**Order Information**

- Specify swing for cupboard doors on **single** units by substituting '\*' with 'R' or 'L' after series number (Units are shown with left swing).
- For pullout shelf, add suffix "-PO" to end of model number (7019L24VAC-PO)

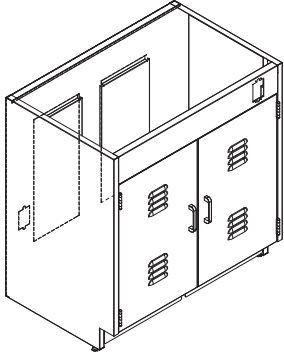

<b>7019*18VAC</b> .....	18"W	<b>7119-35VAC</b> .....	35"W
<b>7019*21VAC</b> .....	21"W	<b>7119-41VAC</b> .....	41"W
<b>7019*24VAC</b> .....	24"W	<b>7119-47VAC</b> .....	47"W
<b>7119-29VAC</b> .....	29"W		



**Narcotic Locker Cabinet**

Model 7601-24NC narcotic locker features vermin-free double wall construction, two adjustable shelves, dual locking double door assembly, and lighted alarm for maximum security of floor stock class narcotic medications. Unit can be wall or counter mounted.

The dual locking double door assembly consists of an 18 gauge double pan welded door and a full height 1/8" plate inside door suspended on a piano hinge welded to the cabinet to eliminate removal of the hinge pins. Both doors are designed with inter-changeable removable core locks\* and are complimented with electronic circuitry. A warning light will automatically be activated if the inside door is left unlocked (the light can be remotely wired to a Director's office to alert entry to the locker). An optional audible alarm is available upon request.

Model **7601-24NC** = 24"W x 30"H x 13"D (standard)

**Modifications for 7601-24NC:**

- Specify swing for cupboard doors on **single** units by substituting '\*' with 'R' or 'L' after series number (Units are shown with left swing).
- Core locks have 50 core locking combinations. Master keying is available.

Narcotics locker wall unit with locking inner & outer doors, lighted alarm assembly.

**7501\*24NC**.....30"H    24"W